

CPAP Mask Catalog



The pause of flow of air intake while sleeping is called as "sleep apnea". This sleeping disorder makes the patient miss one or two breaths while he is asleep. This cycle continues frequently, making the person wake up completely and thereby disrupting the sleep pattern.

Understanding the working principle of the CPAP masks will help you to choose the best comfortable machine for overcoming the sleeping problem. Wearing the CPAP masks during sleep can make people feel uncomfortable, especially if the mask is bulky, hard, heavy, and noisy or causes skin irritation. People using the CPAP masks face five main problems such as air leakage, pressure from headgear, limited sleep positions, skin irritation and noise. However, thankfully, with HSINER's advancements in CPAP mask, all these difficulties have been overcome.

 **HSINER®**
Touching lives, Improving life.



HSINER Co., Ltd.

No. 312, Jhongshan Rd., Shengang Dist., Taichung City 429, Taiwan.

Tel: 886-4-2515-2480 Fax: 886-4-2515-2482

<http://www.hsiner.com> E-mail: sales@hsiner.com

For Pediatric

- P.01~P.02 - *Cirri-mini* - Full Face Mask & Nasal Mask
- P.03~P.04 - **Type 852** - Nasal Mask

For Infant

- P.05 - *NeoQ* - Infant Nasal Mask

For Adult & Pediatric

- P.06~P.08 - **Type 837** - Total Face Mask

For Adult

- P.09~P.10 - **Type 826** - Full Face Mask & Nasal Mask
- P.11~P.12 - **Type 828** - Full Face Mask & Nasal Mask
- P.13~P.14 - *DualSpa™* - Nasal Mask


- P.15~P.16 - **Type 843** - Full Face Mask & Nasal Mask
- P.17~P.18 - **Type 836** - Full Face Mask & Nasal Mask
- P.19~P.20 - **Type 827** - Full Face Mask & Nasal Mask
- P.21~P.22 - **Type 831** - Full Face Mask
- P.23~P.24 - **Type 825** - Full Face Mask & Nasal Mask
- P.25~P.26 - **Type 820** - Full Face Mask & Nasal Mask
- P.27 - **Type 850** - Full Face Mask
- P.28 - **Type 851** - Full Face Mask & Nasal Mask

Pillow Type



- P.29 - *Rolin™ P* - Pillow Nasal Mask
- P.30 - **Type 849** - Pillow Nasal Mask
- P.31 - **Type 101** - Pillow Nasal Mask

FinTek

- FinTek Mask adjustment combines proven performance coupled with the latest nasal and full face mask technology ensuring patient acceptance.
- FinTek enables effective sealing and comfort to enhance patient therapy values and continued patient commitment to treatment.
- FinTek is Available on our full range Nasal and full-face masks maximizing use in a larger cohort of patients, giving confidence to both clinician and patient.

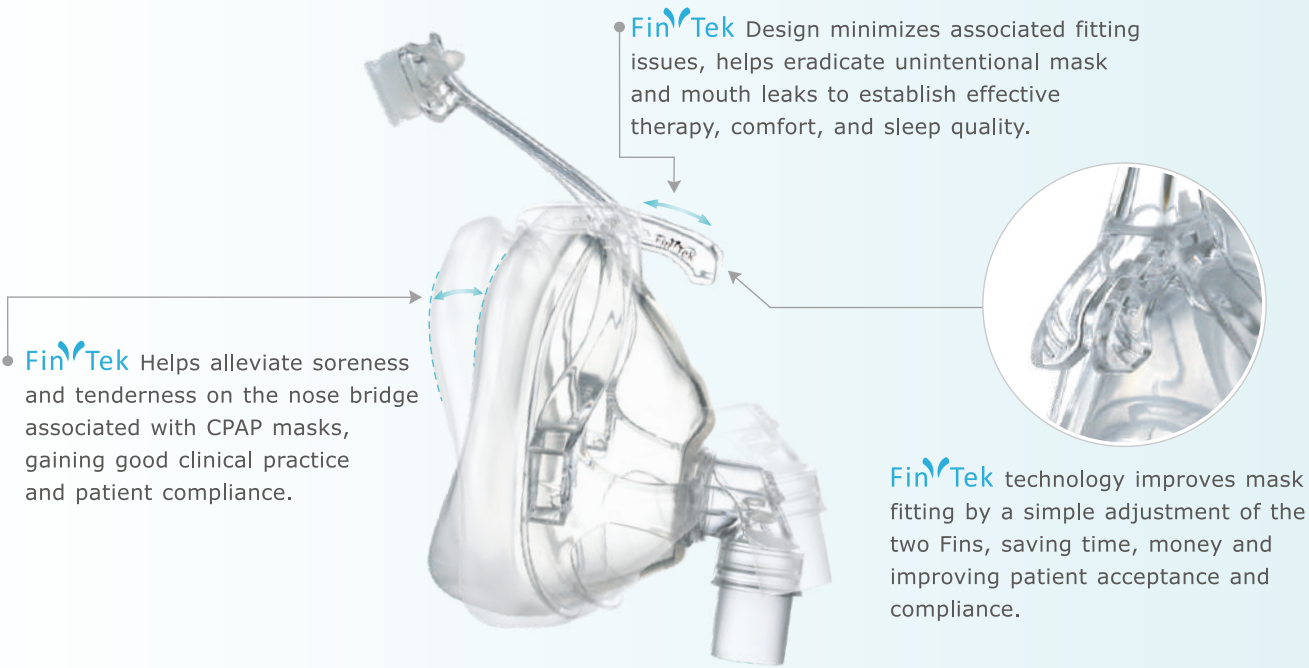


Type *Cirri-mini* Pediatric Full Face & Nasal Mask

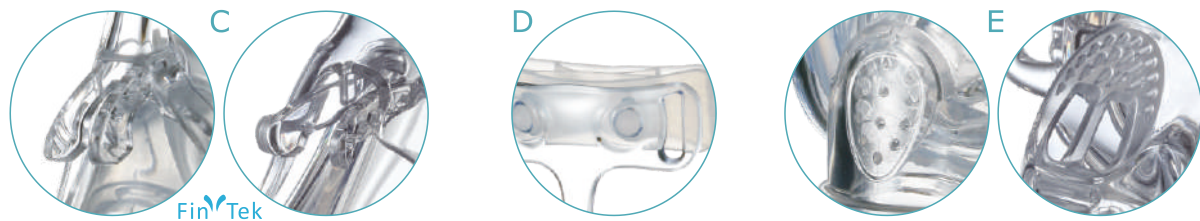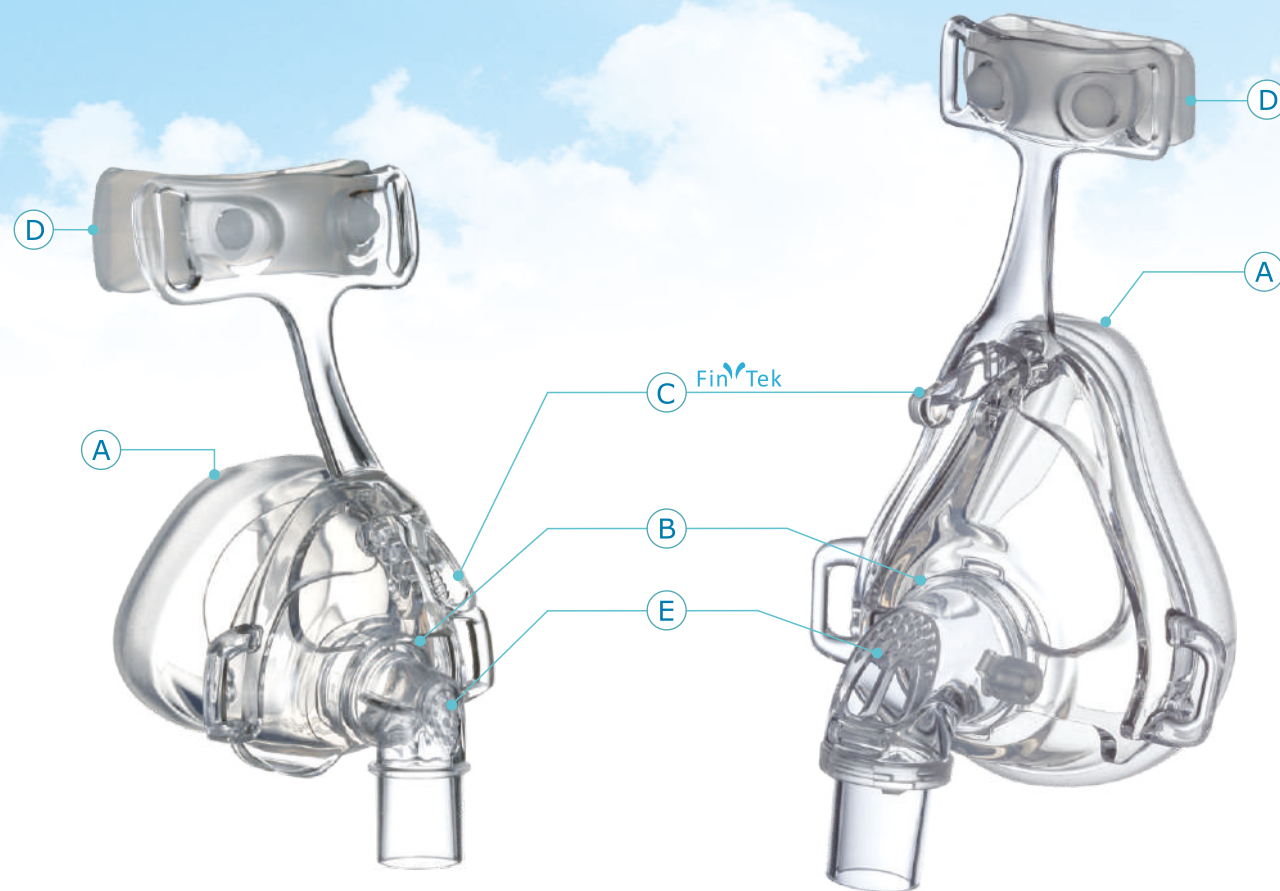

- (A) New 3D mask cushion provides an excellent seal without tight headgear; the face pressure sores and discomfort won't happen.
- (B) Just a simple action to attach and detach the cushion from frame.
- (C) **FinTek** >> Innovative adjustable forehead support provides stability with easy one-handed adjustment and allows the mask fit most of patients.
- (D) Flexible, one-piece forehead pad evenly conforms to the forehead contour, so the pad stays in place and you won't feel a thing.
- (E) Quiet air vent minimizes noise and enable efficient CO₂ washout. Vented air is spread as cloud so that it does not disturb patient or companion.



#11339 : Child S

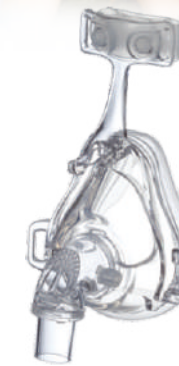

#11340 : Child M



#11341 : Child L



#11342 : Child XL



#11349 : Child S



#11350 : Child M




#11351 : Child L



#11352 : Child XL

Type 852 Pediatric Nasal Mask

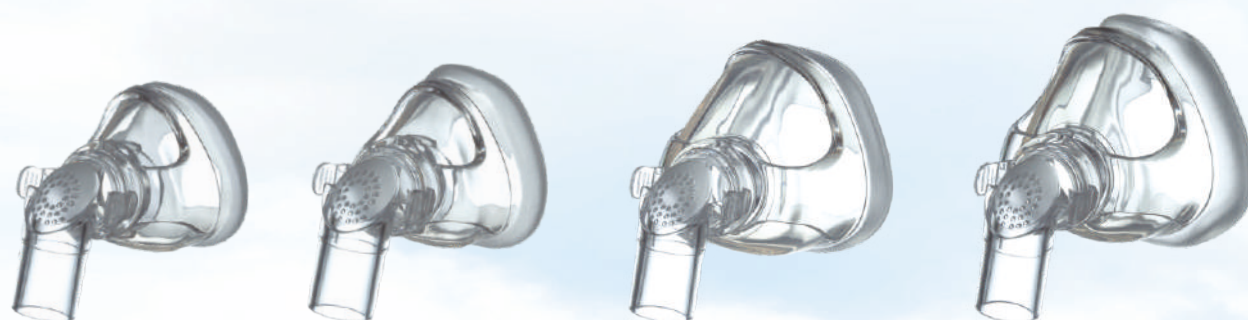

#15349 Child Small
#15359 Child Small


#15350 Child Medium
#15360 Child Medium

#15351 Child Large
#15361 Child Large

#15352 Child X-Large
#15362 Child X-Large

Optional Headgear Color

- Tiffany Blue
- Orange

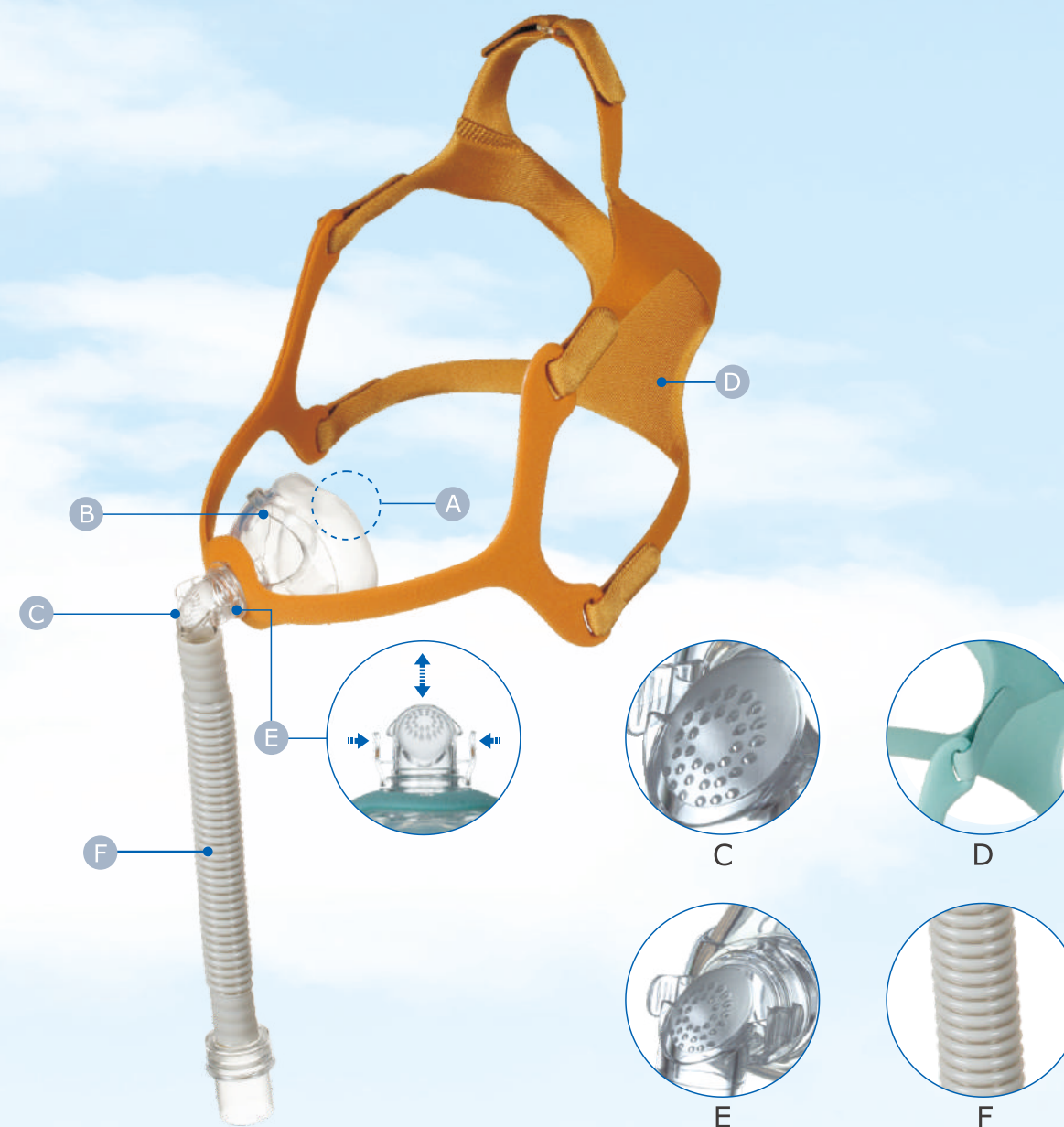

- A New 3D mask cushion provides an excellent seal without tight headgear; the face pressure sores and discomfort won't happen.
- B Just a simple action to attach and detach the cushion from frame.
- C Quite air vent minimizes noise and enable efficient CO₂ washout. Vented air is spread as cloud so that it dose not disturb patient or companion.
- D Optional headgear color available.
- E Quick-Release adapter: Just press the clips for a fast and easy disconnection from the CPAP hose when needed.
- F Flexible Hytrel tube, light weight.

Type

NeoQ

Infant Nasal Mask

Developed to care and support his/her first breath.
The smallest patients often need the most care.
Safety and comfort are the most important to help babies with their earliest steps in life.

NeoQ



Entirely made of liquid soft silicone makes the cushion part perfectly fits the facial characteristics of the patient and ensures minimal pressure points.

The system is designed for baby from 3 to 12 kgs. It consists of a nasal mask, a 5-point headgear and a collapsible tubing.

The silicone mask is available in 4 sizes, and headgear in 2 sizes.



Item No.	Size	Weight Range
13451	#1	<5kg
13452	#2	4~6kg
13453	#3	6~8kg
13454	#4	>8kg

Total Face Mask

Type

837



Total Face Mask
(For Pediatric)
#13830 XX Small



Total Face Mask
(For Pediatric)
#13831 X Small
#13832 Small



Total Face Mask
(For Adult)
#13843 Small - S
#13844 Medium - M
#13845 Large - L
#13846 Extra Large - XL

Type 837 Total Face Mask



- A. The transparent cover creates a crystal clear view for the user.
- B. Silicone total full face cushion fits just right and firmly on the whole face, has less pressure for the user.
 - B-1 : Strengthened stiffness of the cushion decreases leakage, without tightening the headgear, easing the pressure to the user.
 - B-2 : Flexible mask frame fits different facial contours, providing better sealing and comfort.
 - B-3 : Ergonomic mask curve allows user to open mouth with the mask on.
- C. Anti-asphyxia valve, automatically open when the generator pressure fails, enabling patients to breathe by themselves.
- D. Quick release adaptor, press make user release / attach the tubing easily during temporary leaving.
- E. Sampling port allows measuring pressure or ETCO₂.
- F. With nasogastric (NG) port.
- G. Quick release buckles, easy for user to release and take off the mask after sleep.
- H. 360° rotational elbow providing a comfortable tube position.
- I. Applicable to the users who have facial hair and general facial irregularity.



Optional Elbow Connectors



Type **826** Full Face Mask
Nasal Mask




Full Face Mask

#12743 Small - S
#12744 Medium - M
#12745 Large - L

Nasal Mask

#12753 Small - S
#12754 Medium - M
#12755 Large - L

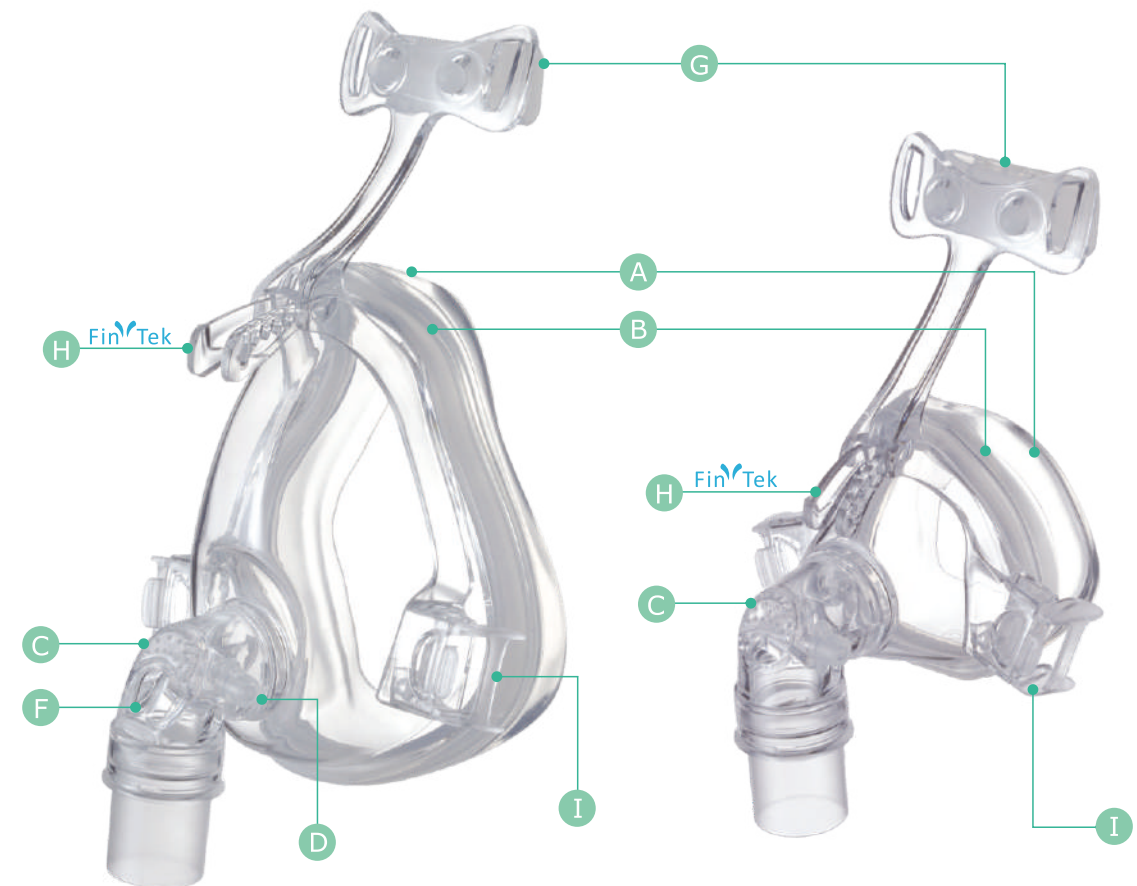

A, B



C, F



G



H

- A** Medical grade liquid silicone full face cushion fits just right and firmly on your face.
- B** Quick assembling design offer user more convenience when taking off / assembling the cushion.
- C** Silence type vent holes provides lower noise and the maximum CO₂ washout, vented air is spread as cloud so that does not disturb patient or companion.
- D** Sampling port allows measuring pressure or ETCO₂.
- E** Minimum components, less components will lower the risk of losing components and easier for user when cleaning.
- F** Anti-asphyxia valve, automatically open when the generator pressure fails, enabling patients to breathe by themselves. (Full Face Mask)
- G** Flexible, one-piece forehead pad evenly conforms to the forehead contour, so the pad stays in place and you won't feel a thing.
- H** FinTek >> Innovative adjustable forehead support provides stability with easy one-handed adjustment and allows the mask fit most of patients.
- I** Push release / attach buckle, easier for user to take off mask after sleep.

Type **828** Full Face Mask
Nasal Mask



Full Face Mask

#12943 Small - S
#12944 Medium - M
#12945 Large - L



Nasal Mask

#12953 Small - S
#12954 Medium - M
#12955 Large - L



- A** Medical grade liquid silicone full face cushion for greater comfort special 3D design for better sealing user's face but don't need to tighten the headgear too hardly.
- B** Quick assembling design offer user more convenience when taking off / assembling the cushion.
- C** Silence type vent holes provides lower noise and the maximum CO₂ washout, vented air is spread as cloud so that does not disturb patient or companion.
- D** Sampling port allows measuring pressure or ETCO₂.
- E** Minimum components, less components will lower the risk of losing components and easier for user when cleaning.
- F** Anti-asphyxia valve, automatically open when the generator pressure fails, enabling patients to breathe by themselves. (Full Face Mask)
- G** Flexible, one-piece forehead pad evenly conforms to the forehead contour, so the pad stays in place and you won't feel a thing.
- H** FinTek >> Innovative adjustable forehead support provides stability with easy one-handed adjustment and allows the mask fit most of patients.
- I** Push release / attach buckle, easier for user to take off mask after sleep.
- J** 360° rotational ball in elbow.

Type *DualSpa*TM Nasal Mask




Nasal Mask

#13253 Small - S

#13254 Medium - M

#13255 Large - L

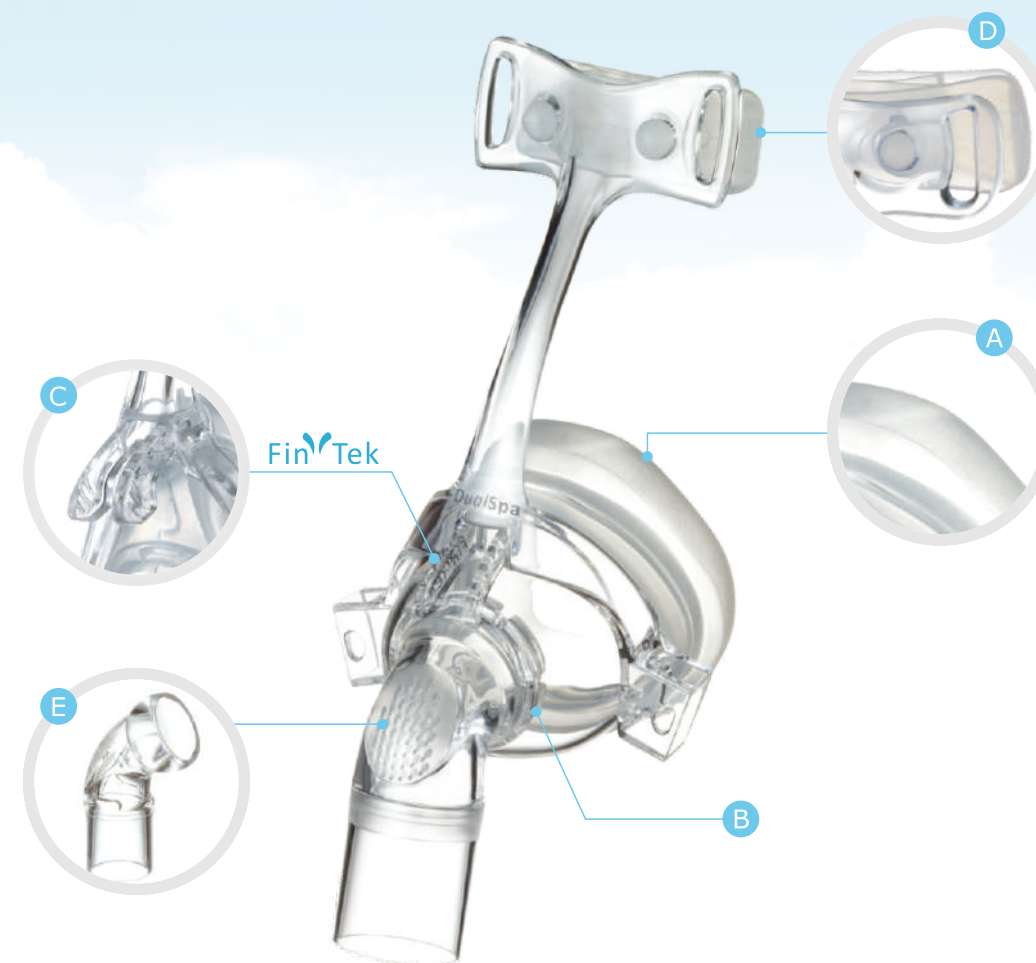

- A** New 3D mask cushion provides an excellent seal without tight headgear; the face pressure sores and discomfort won't happen.
- B** Just a simple action to attach and detach the cushion from frame.
- C** *FinTek* >>Innovative adjustable forehead support provides stability with easy one-handed adjustment and allows the mask fit most of patients.
- D** Flexible, one-piece forehead pad evenly conforms to the forehead contour, so the pad stays in place and you won't feel a thing.
- E** 360°rotational ball in elbow.

Type

843

Full Face Mask

Nasal Mask




Full Face Mask

- #14443 Small - S
- #14444 Medium - M
- #14445 Large - L

Nasal Mask



- #14453 Small - S
- #14454 Medium - M
- #14455 Large - L

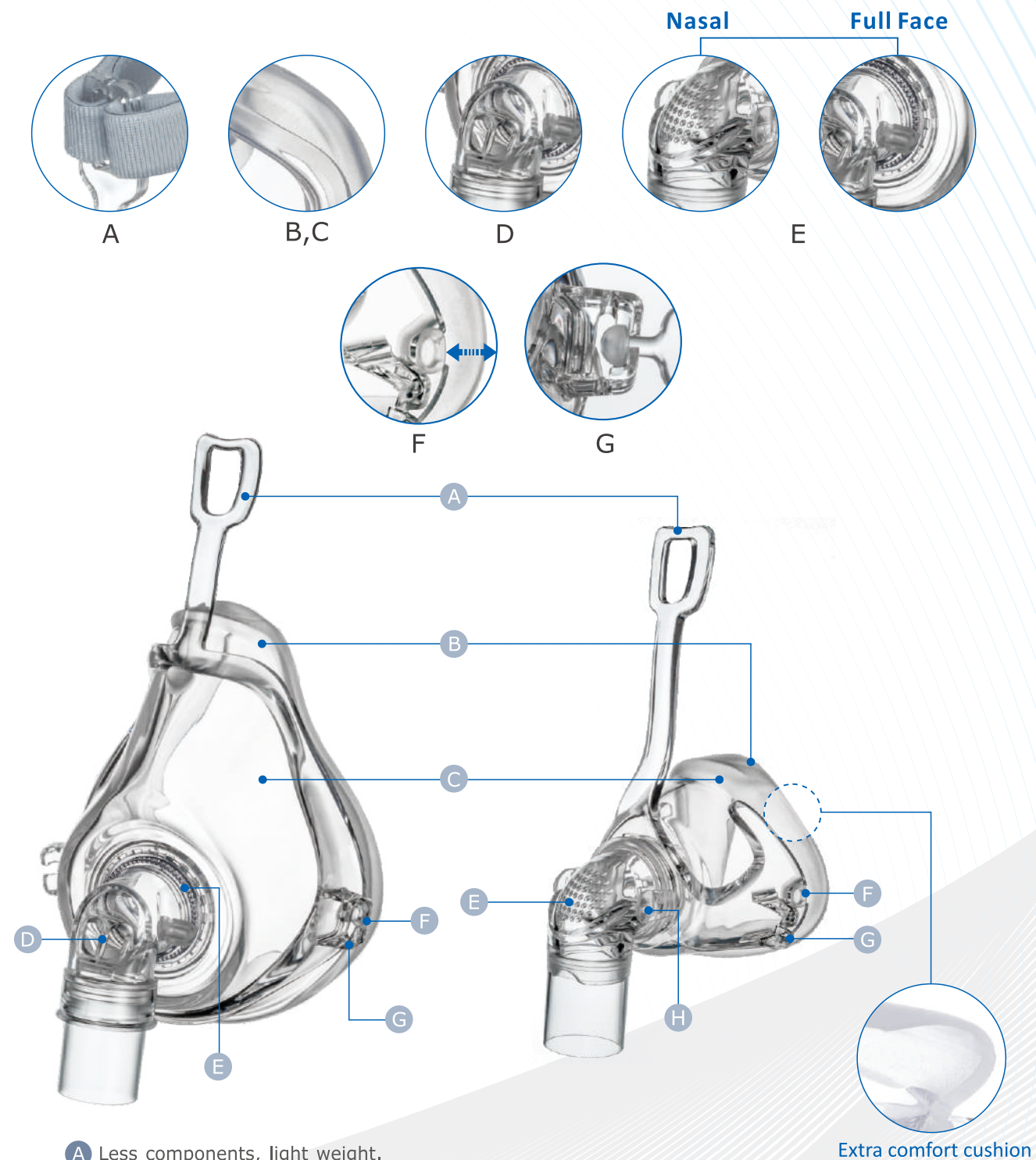

- A Less components, light weight.
- B Anatomically-contoured Comfort for fit and seal.
- C Medical liquid silicone construction.
- D Anti-asphyxia valve, automatically open when the generator pressure fails, enabling patients to breathe by themselves.
- E Silence type vent holes provides lower noise and the maximum CO2 washout, vented air is spread as cloud so that does not disturb patient or companion.
- F Latch on design for easy assembling: using two protruding knob on either side to latch on to the gap on the frame.
- G Ball-type design for easy assembling/disassembling.
- H Quick release elbow connector, easy to take off if users need to have sudden depart from the device.

Extra comfort cushion

Type **836** Full Face Mask
Nasal Mask

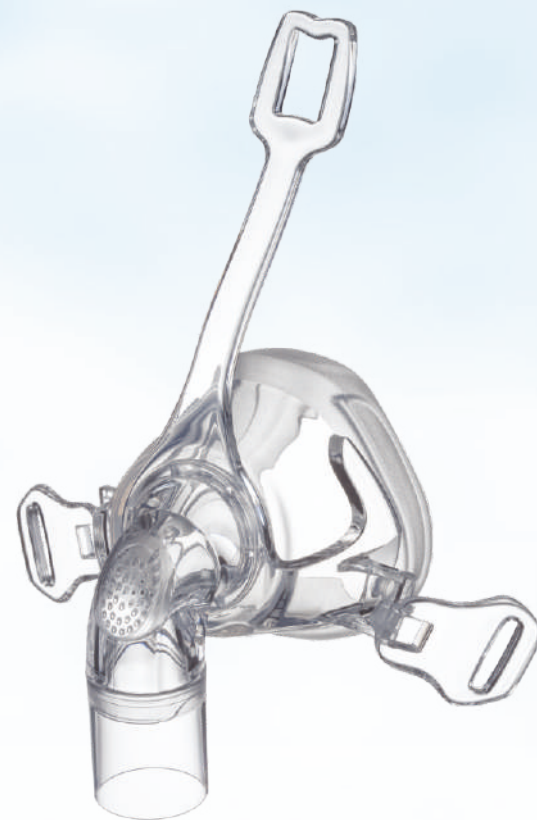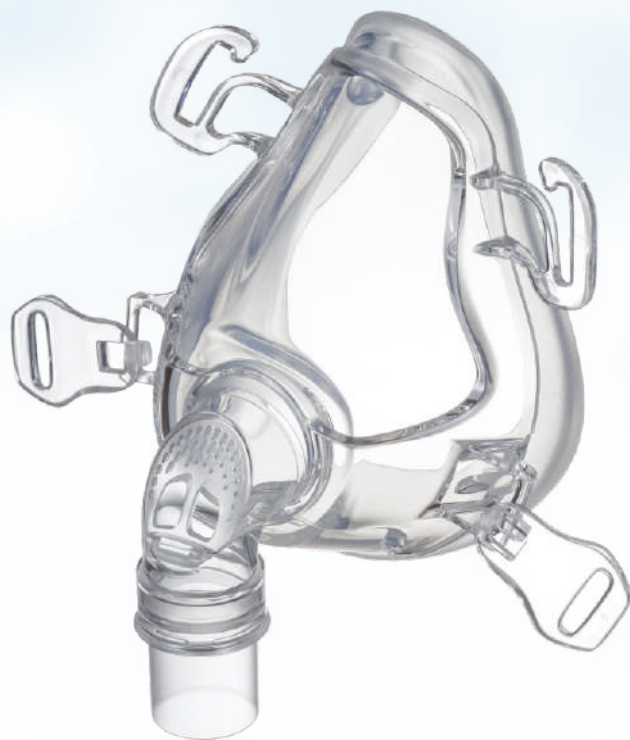


Full Face Mask

#13742 Extra Small - XS
#13743 Small - S
#13744 Medium - M
#13745 Large - L



Nasal Mask


#13753 Small - S
#13754 Medium - M
#13755 Large - L

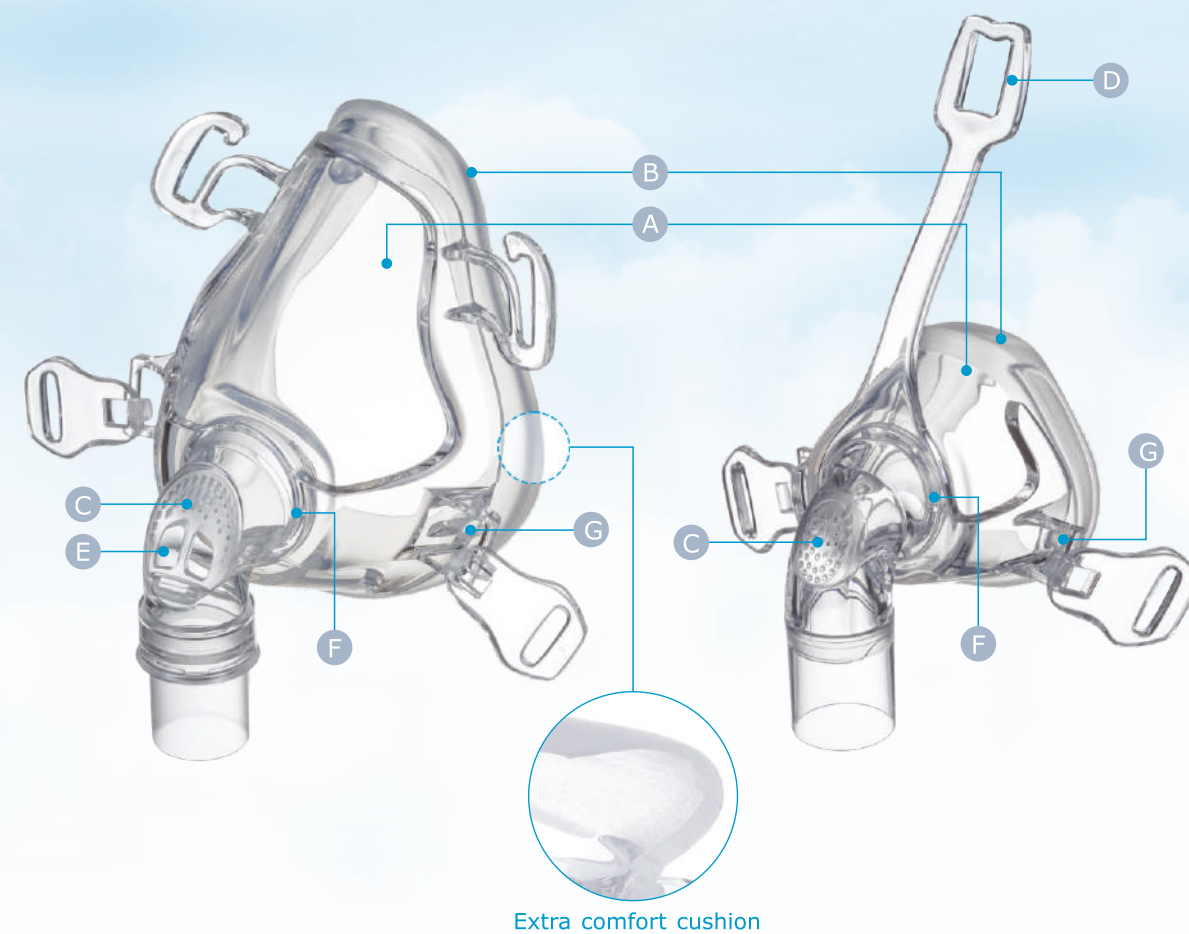


A, B



C, E

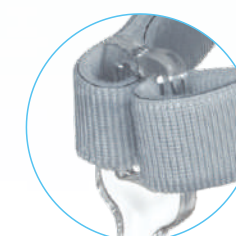

D

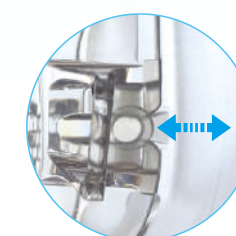


G

- A Medical liquid silicone construction.
- B Anatomically-contoured Comfort for fit and seal.
- C Silence type vent holes provides lower noise and the maximum CO₂ washout, vented air is spread as cloud so that does not disturb patient or companion.
- D Less components, light weight.
- E Anti-asphyxia valve, automatically open when the generator pressure fails, enabling patients to breathe by themselves.
- F 360 degree + floating swivel L connector provides greater flexibility and less disturbance while rolling over in sleep, in order to further advance your sleep quality.
- G Latch on design for easy assembling: using two protruding knob on either side to latch on to the gap on the frame.

Type 827 Nasal Mask



Item No.
#12853 Small/Medium
#12855 Large

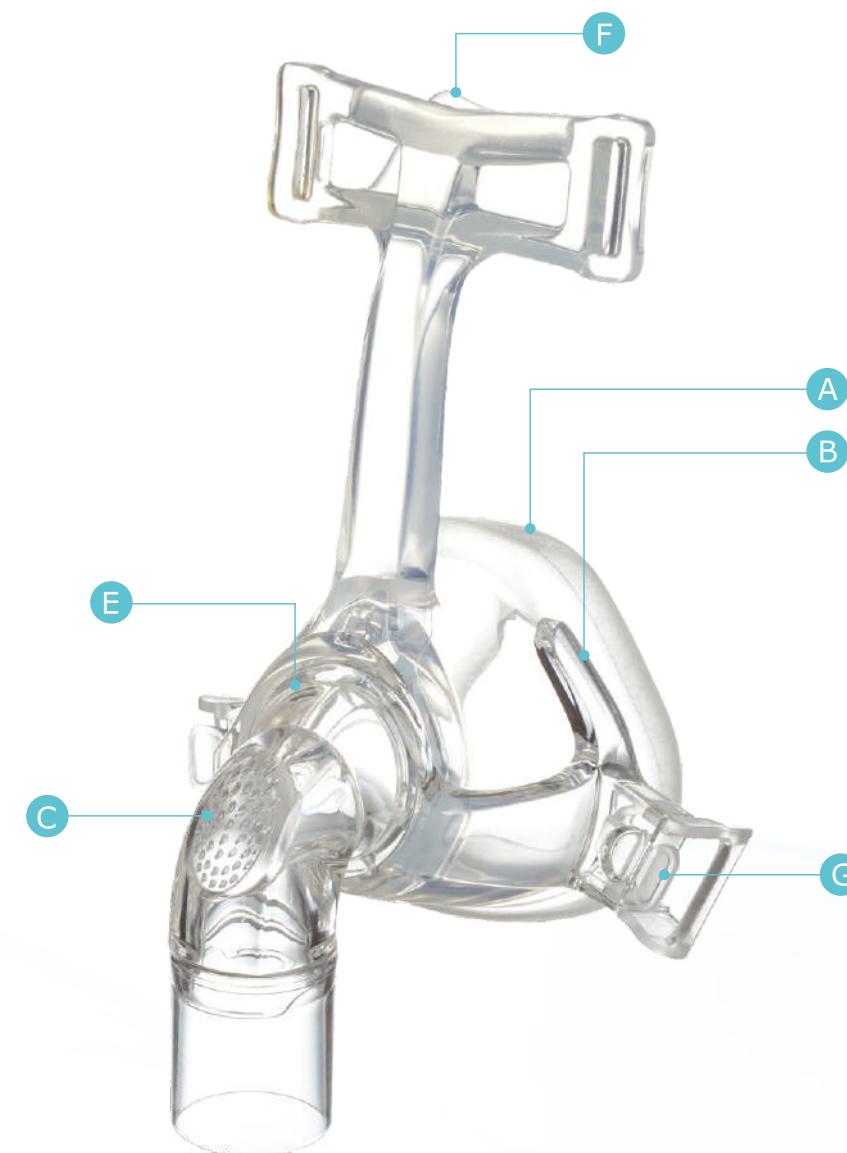

A, B



C, E



F




G

- A** Integrate medical grade liquid silicone cushion provides the greatest comfort to users.
- B** Quick assembling design offer user more convenience when taking off / assembling the cushion.
- C** Silence type vent holes provides lower noise and the maximum CO₂ washout, vented air is spread as cloud so that does not disturb patient or companion.
- D** Minimum components, less components will lower the risk of losing components and easier for user when cleaning.
- E** 360 degree + floating swivel L connector provides greater flexibility and less disturb as turn over while sleeping in order to advance further your sleep quality.
- F** Flexible, one-piece forehead pad evenly conforms to the forehead contour, so the pad stays in place and you won't feel a thing.
- G** Push release / attach buckle, easier for user to take off mask after sleep



Full Face Mask

- #13243 Small - S
- #13244 Medium - M
- #13245 Large - L

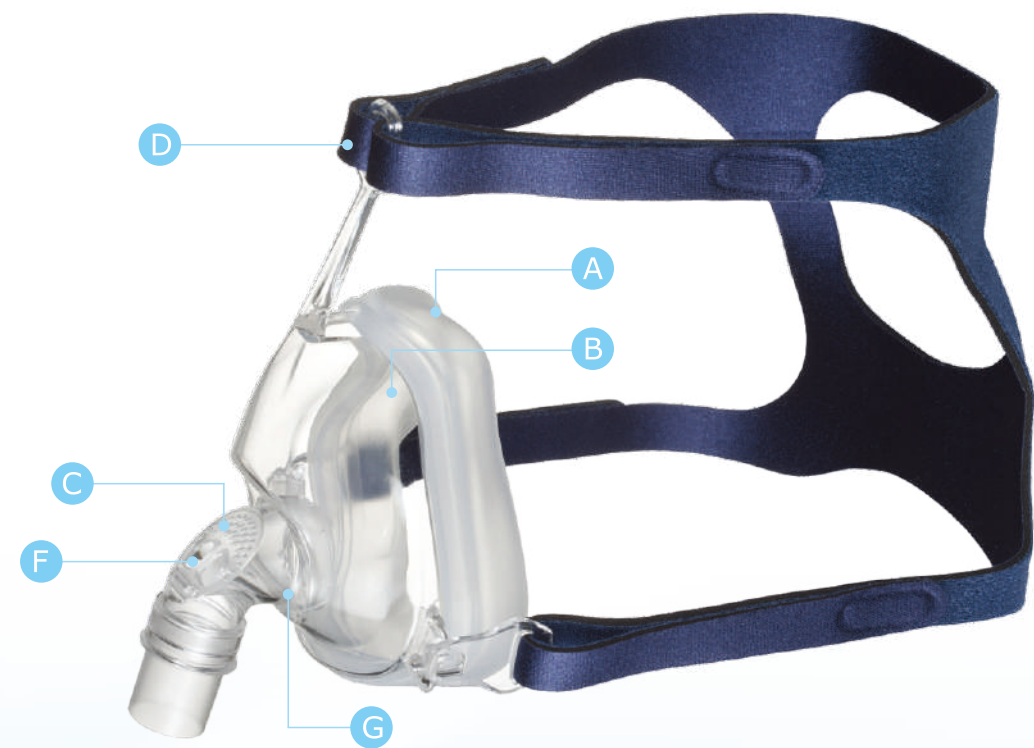

- A** Medical grade liquid silicone full face cushion for greater comfort special 3D design for better sealing user's face but don't need to tighten the headgear too hardly.
- B** Quick assembling design offer user more convenience when taking off / assembling the cushion.
- C** Silence type vent holes provides lower noise and the maximum CO₂ washout, vented air is spread as cloud so that does not disturb patient or companion.
- D** Less components, light weight.
- E** Minimum components, less components will lower the risk of losing components and easier for user when cleaning.
- F** Anti-asphyxia valve, automatically open when the generator pressure fails, enabling patients to breathe by themselves.
- G** 360° rotational ball in elbow.





Full Face Mask

#12643 Small - S
#12644 Medium - M
#12645 Large - L



Nasal Mask

#12653 Small - S
#12654 Medium - M
#12655 Large - L

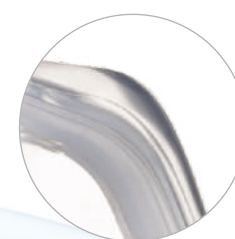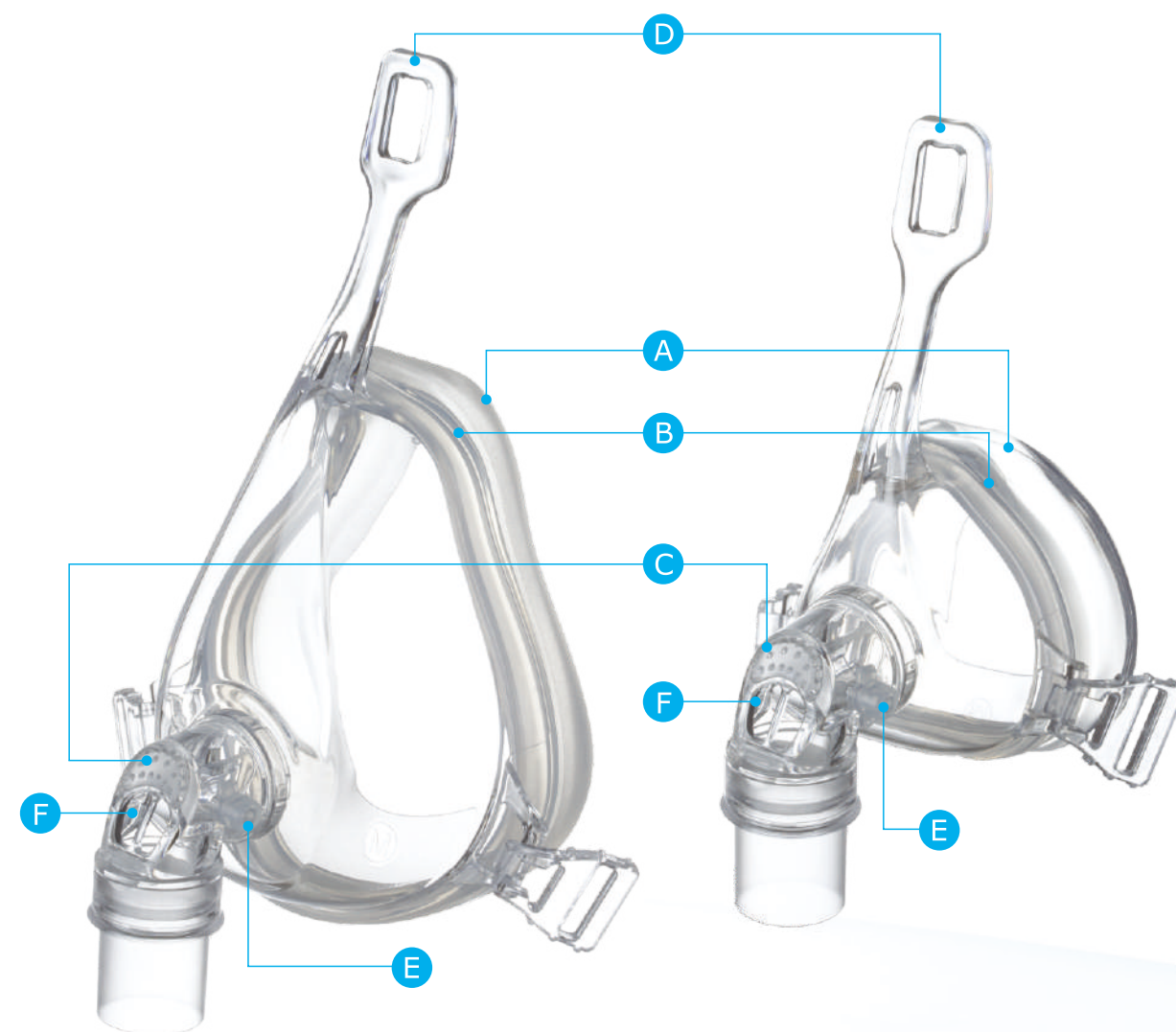

A, B



C, F



D

- A Medical grade liquid silicone full face cushion fits just right and firmly on your face.
- B Quick assembling design offers users more convenience when taking off / assembling the cushion.
- C Silence type vent holes maximize CO₂ washout and reduce noises. Also, the vented air is spread like a cloud, creating fewer disturbances for patient or companion.
- D Less components, light weight.
- E Sampling port for pressure and ETCO₂ measurement at the same time.
- F Anti-asphyxia valve, automatically opens when the generator pressure fails, enabling patients to breathe by themselves. (Full Face Mask)
- ✂ Minimum components allow lower risk of losing components and easier for cleaning.





Full Face Mask

#12143 Small - S
#12144 Medium - M
#12145 Large - L

Nasal Mask

#12153 Small - S
#12154 Medium - M
#12155 Large - L

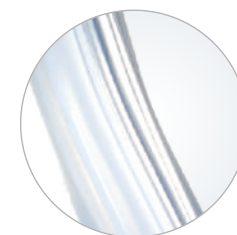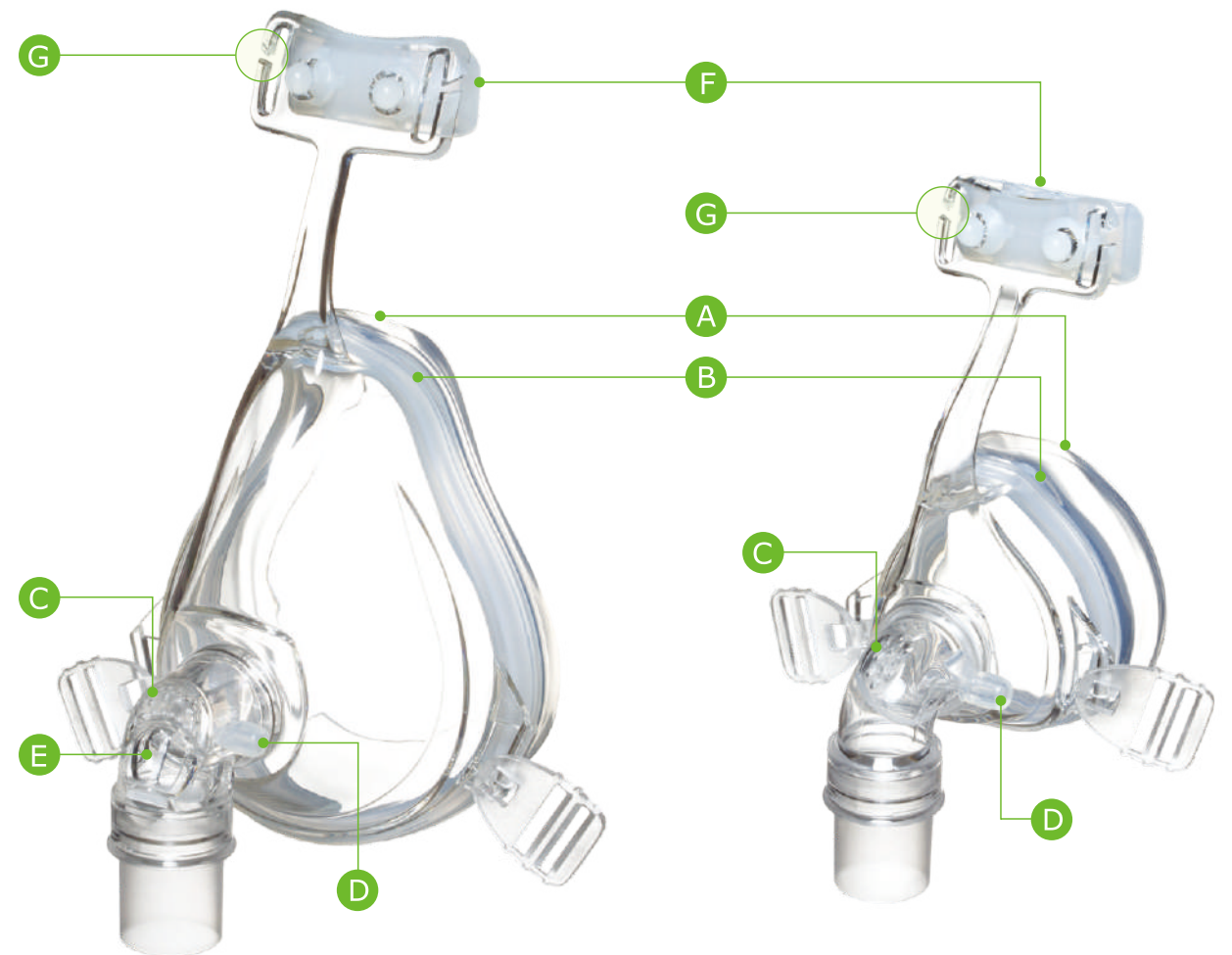

A,B



C,E




F,G

- A** Medical grade liquid silicone full face cushion for greater comfort special 3D design for better sealing user's face but don't need to tighten the headgear too hardly.
- B** Quick assembling design offer user more convenience when taking off / assembling the cushion.
- C** Silence type vent holes provides lower noise and the maximum CO₂ washout, vented air is spread as cloud so that does not disturb patient or companion.
- D** Sampling port allows measuring pressure or ETCO₂.
- E** Anti-asphyxia valve, automatically open when the generator pressure fails, enabling patients to breathe by themselves.(Full Face Mask)
- F** Flexible, one-piece forehead pad evenly conforms to the forehead contour, so the pad stays in place and you won't feel a thing.
- G** Easy to unhook or reattach the upper headgear by threading it through the loop.



Type **Rolin™ P** Nasal Pillow Mask

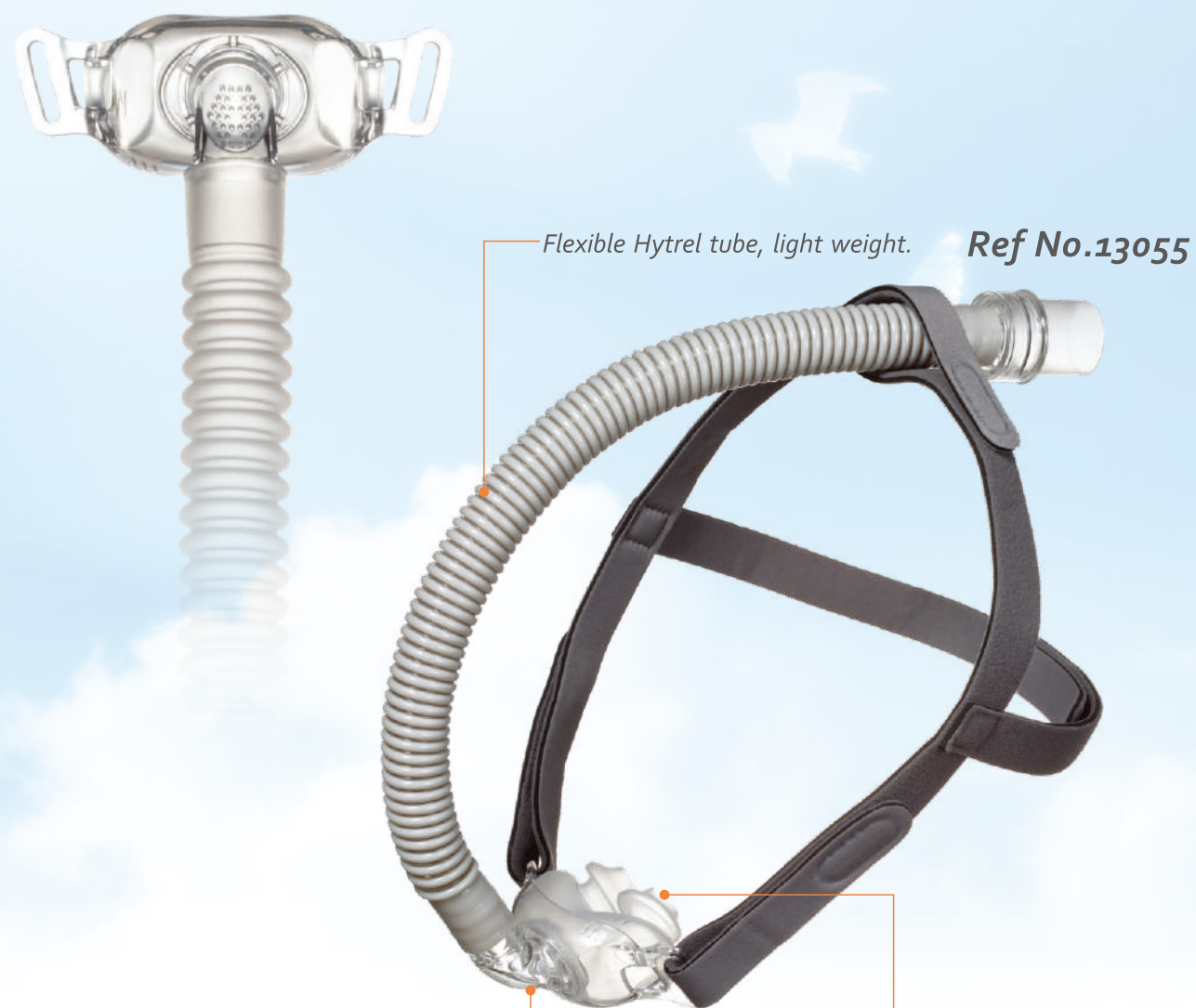

360° rotational ball in elbow.




3 different pillow sizes, one piece design, easy to assemble and disassemble.



Nasal Pillow Mask Type **849**

A simple comfort that you never experienced.

Ref No.15055

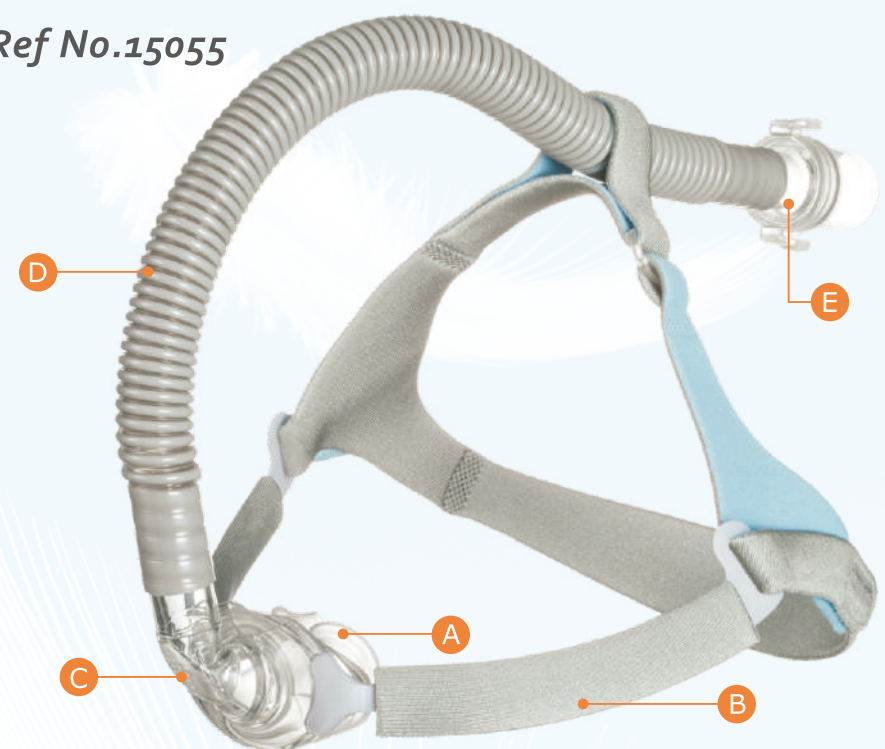

Size: L

Size: M

Size: S



A

E

- A** 3 different pillow sizes, one piece design, easy to assemble and disassemble.
- B** Soft sheaths which are comfortable and fit to the skin without pressure mark.
- C** 360° rotational ball in elbow.
- D** Flexible Hytrel tube, light weight.
- E** Quick-Release adapter: Just press the clips for a fast and easy disconnection from the CPAP hose when needed.

Type 101 Nasal Pillow Mask

CAPA Nasal pillow Mask have been well received by patients who want something innovative and different.

In general a nasal pillow uses two soft silicone anatomically shaped "Pillows" which are inserted slightly into the nostril. Because of the unique designs of HSINER such as 4 different pillow sizes and one-piece injected silicone frame, the masks can be quite comfortable.

There is minimal facial contact and allow for almost any sleeping position without mask displacement. The minimal and lightweight design is welcomed by CAPA mask users over the traditional bulky CPAP nasal or full face mask.

Additionally, the exceptionally quiet operation of the nasal pillow helps the efficient CO₂ washout and the sleep experience even at higher pressures.

Ref No.11254



4 different sizes prong
& 4 different pitches

Features:

- ▶ Easy to fit, assemble, and maintain.
- ▶ 4 different pillow sizes(XS. S. M. L) provide diverse compliance of the pitch of two nostrils.
- ▶ 360° rotational elbow with flexible silicone tube.
- ▶ Quiet vent design.